




Relaxing music
Built-in Bluetooth music function



Comfortable chair and relaxing music

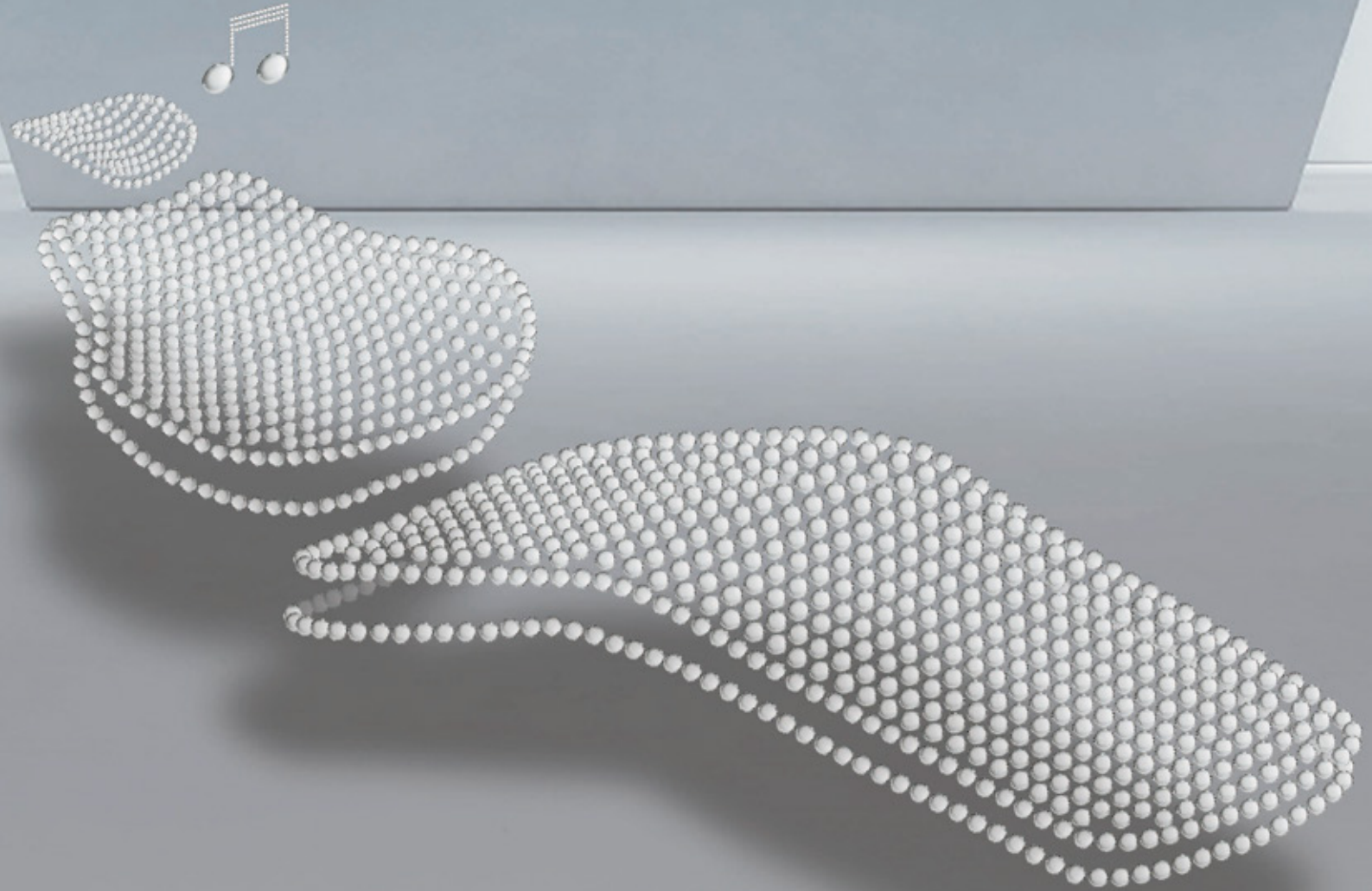

Hanging Type



Top mounted



Cart Type





Bluetooth music function
Relaxing music



Ceramic spittoon
180 degrees rotatable, deep and large capacity bowl.



Rotatable assistant tray by 180 degrees



European style swing type dentist tray

Silicon handpiece hose pad and instrument tray
Flexible and easy to use



Hanging type dentist tray

Acrylic protective cover
Sterilisable silicone handrail cover

Cingol 3.0 Smart Operation System

- Self-diagnostic detection system.
- Emergency stop and first-aid chair position.
- Interlock system: when the handpiece are being used, the chair will not move.




Touch panel dentist tray



Ti-Alloy handpiece connector



Imported 3-way syringe

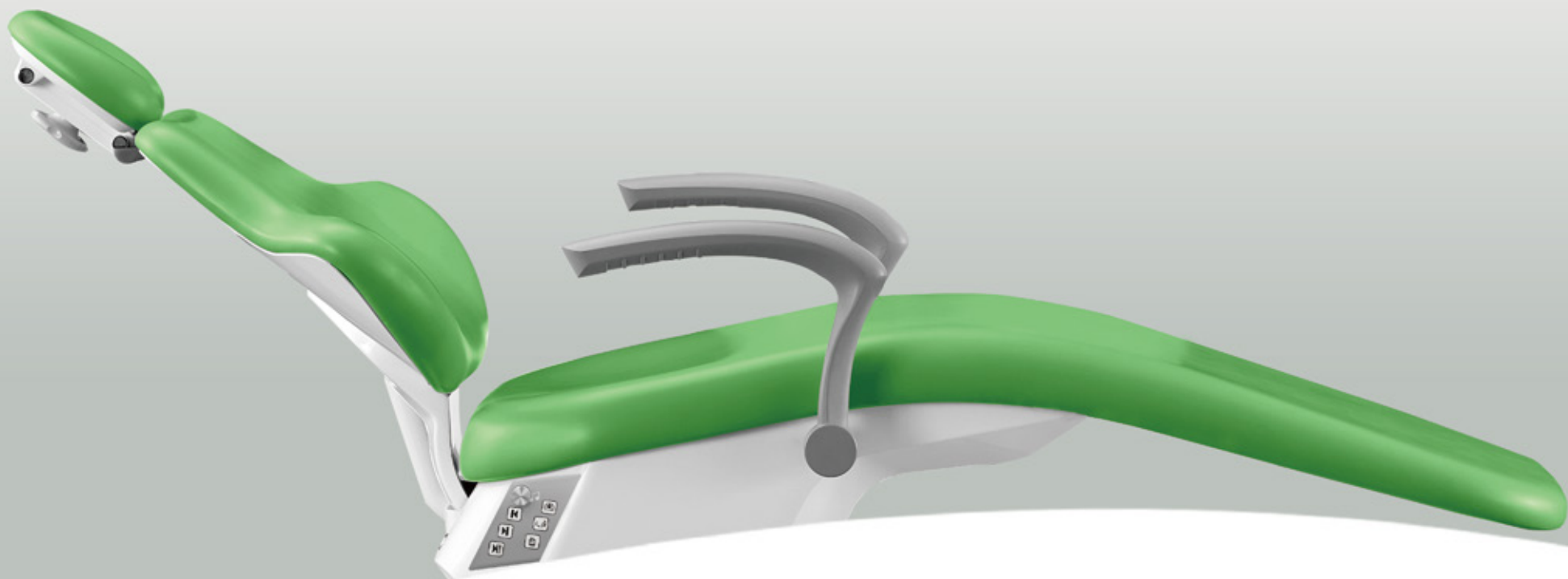


Ergonomic Upholstery

The thin backrest allows more room for the dentist.



Slip resistant armrest

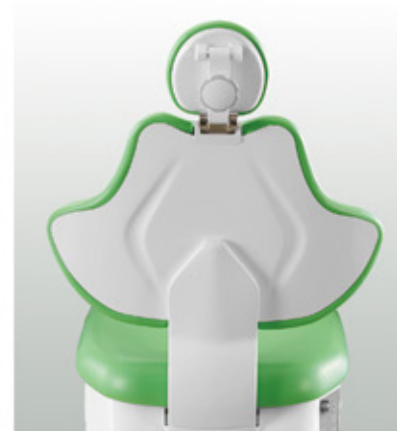

Metal backrest

Detachable backrest for easy replacement of cushion.



3-axis headrest

Can be positioned for adults and children or disabled alike.



70 degrees rotatable sidebox



Detachable sidebox cover enables easy maintenance.




Rotary plug-in double sealed filter, external design, easy to clean.



Built in double bottle water storage system for distilled water. The water pours straight into the bottles from the top of the sidebox.

T.MASTER

Stepless adjustment of shadowless cold light,
close to natural light.

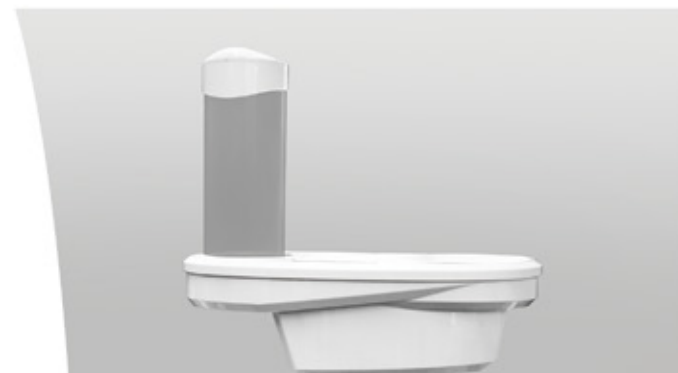

3 in 1 tissue holder with universal
USB charging port.



The multi-function pedal also
controls the light and cup filling.



The Italian style light arm is light
weight and flexible.

OPTIONS



LED curing light



Scaler



Handpiece



Electric motor



Intra oral camera with support



ECG



T-Master dental simulator




T-Master teaching recording lamp



WINQ oilless air compressor



WinQ wet suction system

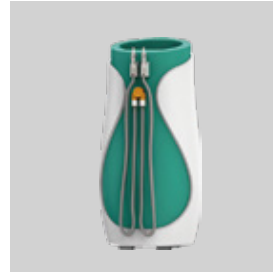


WinQ portable air suction




WinQ air suction system

SPECIFICATION


Input voltage: 230V 10%,50Hz
 Water pressure: 0.2-0.4Mpa/2.0-4.0 bar
 Air pressure: 0.55-0.8Mpa/5.5-8.0 bar
 Gross Weight: 275kg

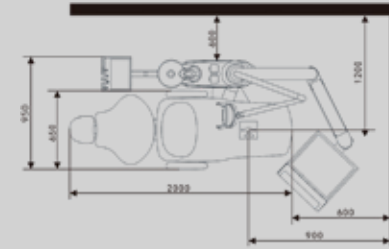

Built-in floor-box installation sketch map

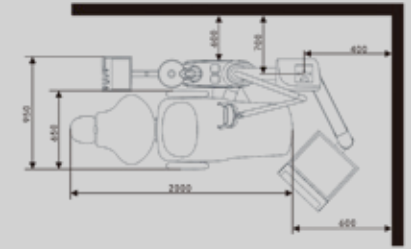

External floor-box installation sketch map

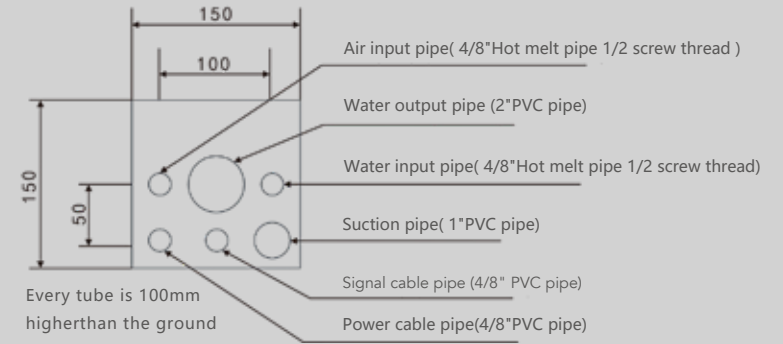

Pre-buried pipe distribution map

Upholstery colours

PU

* Please refer to actual color swatches for upholstery leather colors.



Popular colours

Options: Fibre leather





CINGOL

Concern about International dental Needs,
our Goal is your Obliging Laugh.



Remark:
We reserve the right to improve all the products
without notice in advance;
The product's picture is only for your reference;
the copyright is reserved and belongs to Cingol.

FOSHAN CINGOL MEDICAL INSTRUMENT CO., LTD
Tel:0086-757-81797099 Fax:0086-757-81797055
Email:cingol@cingol.com QQ:2851731555
Site:www.cingol.com
Add: Building 6-1, Nansha Industrial Village, Jiangli North Rd,
Danbao Town, Nanhai District, Foshan, China.