

Modern design  
Get yourself an outstanding dental chair





Hanging Type



Top Mounted



Cart Type





**Bluetooth Music Function**  
Relaxing music



**Detachable spittoon**  
180 degrees rotatable



**Metal backrest**  
Easy cushion replacement



**Disinfection**  
Easy operation of automatic disinfection of the handpiece hoses  
and suction hoses



- European style dentist tray
- Ease of operation



- Cart type unit provides freedom of movement.**
- Adjustable Height
  - Emergency stop and first-aid chair position
  - Clean-up chair position

## Cingol 3.0 Smart Operation System

- Self-diagnostic detection system.
- Emergency stop and first-aid chair position.
- Interlock system: when a handpiece or scaler is being used, chair movement (up, down, back, forward) is locked.



Touch sense control panel



Ti-alloy handpiece tubing connector



Imported-brand 3-way syringe



## Ergonomic Upholstery

The thin backrest allows more room for the dentist



Multi-joint assistant tray arm.  
Multi-angle movement convenient  
for four hands operation.



Touch sense assistant tray



3 Axis can be positioned for adults and  
children or disabled alike





Sidebox 90 degrees rotatable



▶ Detachable sidebox cover, easy maintenance.



▶ Filter in each suction tube. Easy cleaning.





▶ Built-in double bottle water storage system for distilled and disinfected water.

# T.MASTER

European style LED  
operating light

Stepless adjustment of LED light level

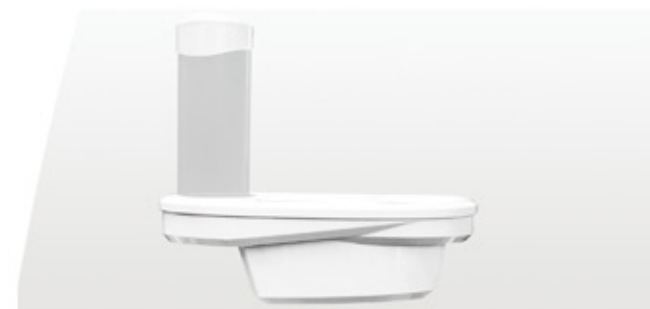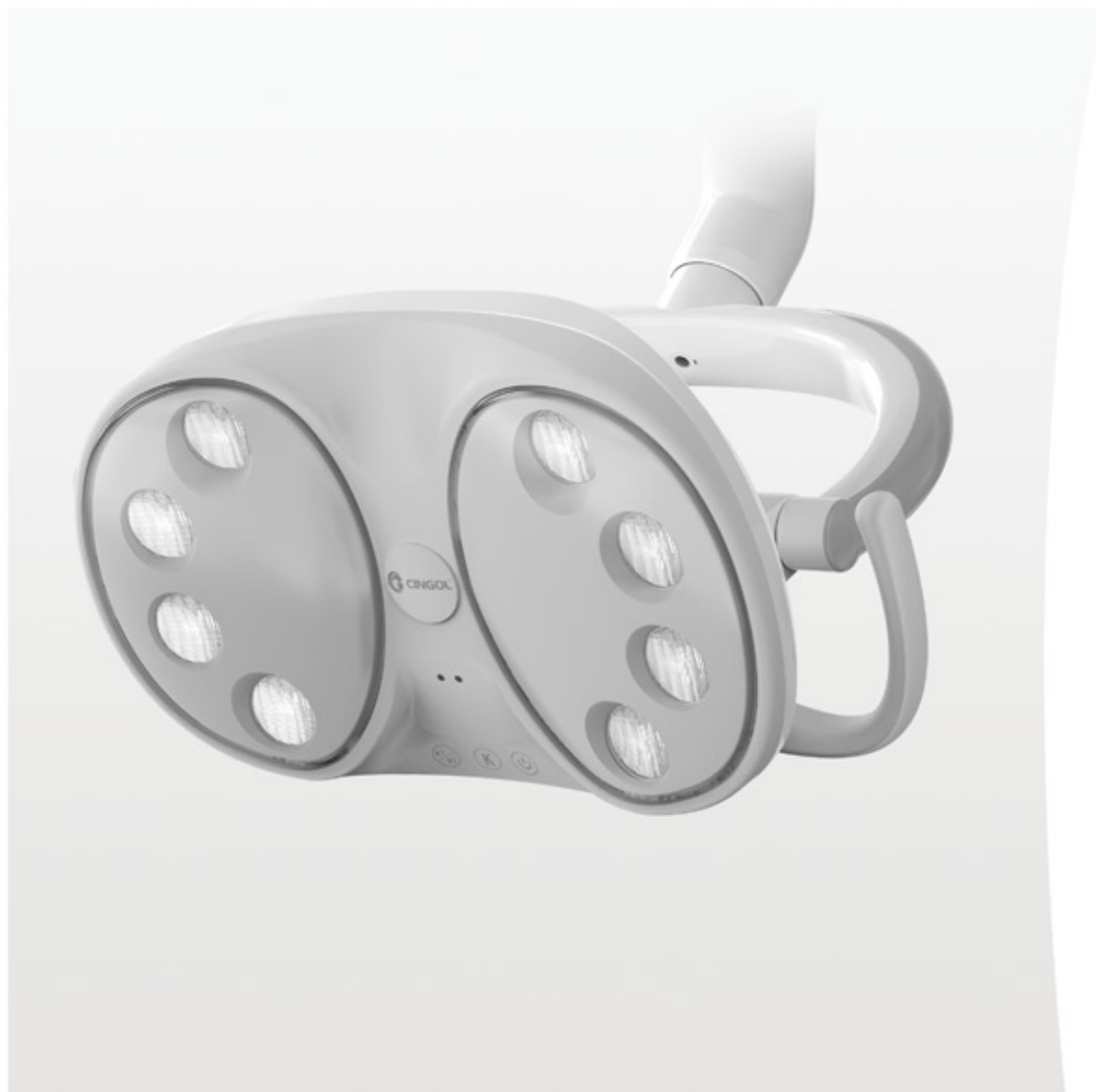

3 in 1 tissue holder with universal  
USB charging port.



Luxury multi-functional foot pedal also controls  
the light and cup filling.



WINQ brand, the Italian style light arm is  
flexible and light.

# OPTIONS



LED curing light



Scaler



Handpiece



Electric motor



Intra oral camera with support



ECG



T-Master dental simulator




T-Master teaching recording lamp



WINQ oilless air compressor



WinQ wet suction system

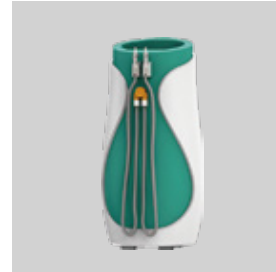


WinQ portable air suction




WinQ air suction system

# SPECIFICATION


Input voltage: 230V 10%,50Hz  
 Water pressure: 0.2-0.4Mpa/2.0-4.0 bar  
 Air pressure: 0.55-0.8Mpa/5.5-8.0 bar  
 Gross Weight: 275kg

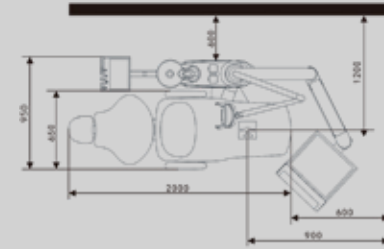

Built-in floor-box installation sketch map

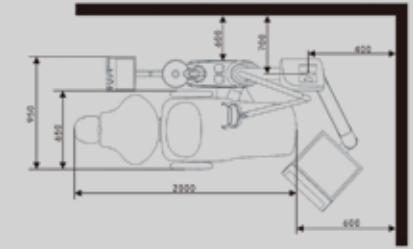

External floor-box installation sketch map

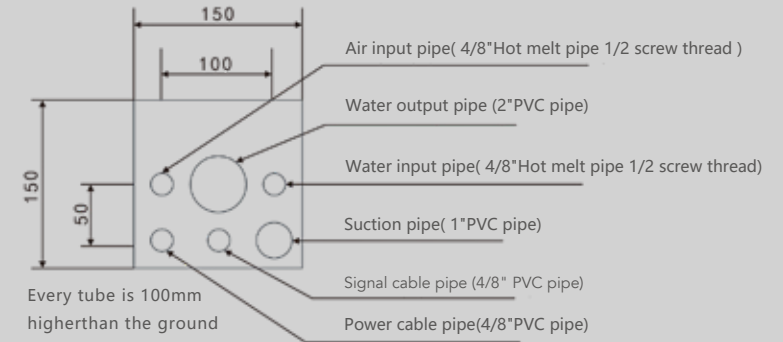

Pre-buried pipe distribution map

## Colour options

\* Please refer to actual color swatches for upholstery leather colors.



> Popular colours ★★★

★★★ S01 ★★ S06 ★ S09

Fibre leather



# CINGOL

Concern about International dental Needs,  
our Goal is your Obliging Laugh.



Remark:

We reserve the right to improve all the products  
without notice in advance;

The product's picture is only for your reference;  
the copyright is reserved and belongs to Cingol.

FOSHAN CINGOL MEDICAL INSTRUMENT CO., LTD  
Tel:0086-757-81797099 Fax:0086-757-81797055  
Email:cingol@cingol.com QQ:2851731555  
Site:www.cingol.com  
Add: Building 6-1, Nansha Industrial Village,Jiangli North Rd,  
Danzao Town,Nanhai District,Foshan,China.