



Extraordinary design  
Focus on craftsmanship



External and internal cultivation is  
real superiority



Hanging Type



Top mounted



Cart Type



# Cingol 3.0 Smart

## Operation System

- Self-diagnostic detection system.
- Shock-resistant function.
- Clean up chair position, one touch to turn off the light, the chair then raises to the highest position and the spittoon will flush automatically.
- Safety emergency stop, first-aid chair position.
- Interlock system, when a handpiece or scaler is being used, chair movement (up, down, back, forward) is locked.



### C-LINK Touch-Screen System



The C-link touch screen system has the built-in functions of electronic micromotor control, memory positions, air pressure and chair movement.



#### Autocup

When a cup is placed under the spout, it automatically fills.



#### Alarm Button

For patients to attract the attention of the dentist.

# Comfortable for the patients.

Bluetooth music function.

The upholstery is fibre leather with high density memory foam.



Excellent Ergonomic Upholstery

# Ergonomically designed upholstery



- Automatic seat positioning for spittoon use.
- Water temperature adjustable for all handpiece tubes.
- The thin backrest gives more room for the dentist to move around and operate.

Metal backrest 70 degrees sidebox rotation  
Detachable backrest for easy cushion replacement.



Touch sense assistant tray  
2 drawers for disposables



Easy to remove and replace



Top-mounted dentist's tray



Ti-alloy handpiece tubing connector



## Elegant design

- 9 memory programs.
- Rotatable 6 rack handpiece holder.
- Digital display of air pressure, stop watch and water temperature.



# Hygiene and Disinfection



- Detachable sidebox cover, easy for maintenance.
- One touch spittoon cleaning and smart disinfection.



Hose disinfection protects dentists and patients from cross infection. Built in double bottle water supply system for distilled and disinfected water. The water is poured into the bottles from the top of the sidebox. Easy to use. Two built-in disinfecting bottles design ensures fast flushing of the handpieces after disinfection.



The automatic ASTER system, which is controlled by one button, provides disinfection for all handpiece hoses and suction hoses.




Easy to access suction line filters. The spittoon has one button disinfection and is detachable.

# T.MASTER

European style LED  
operating light

Stepless adjustment of shadowless cold light, close to natural light.

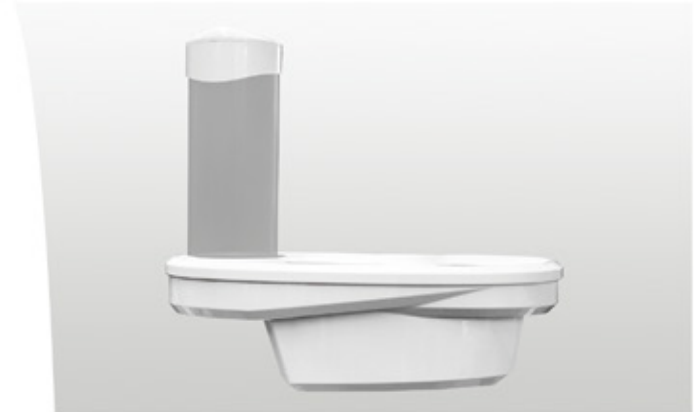

3 in 1 tissue holder with universal  
USB charging port.



The multi-function pedal also controls  
the light and cup filling.



WINQ brand, the Italian style light  
arm is flexible and light.

# OPTIONS



LED curing light



Scaler



Handpiece



Electric motor



Intra oral camera with support



ECG



T-Master dental simulator




T-Master teaching recording lamp



WINQ oilless air compressor



WinQ wet suction system

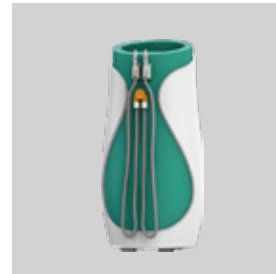


WinQ portable air suction




WinQ air suction system

# SPECIFICATION


Input voltage: 230V 10%,50Hz  
 Water pressure: 0.2-0.4Mpa/2.0-4.0 bar  
 Air pressure: 0.55-0.8Mpa/5.5-8.0 bar  
 Gross Weight: 285kg

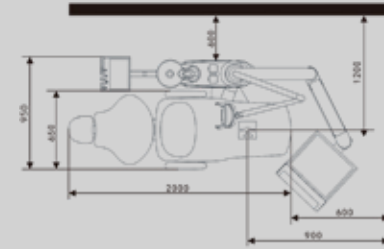

Built-in floor-box installation sketch map

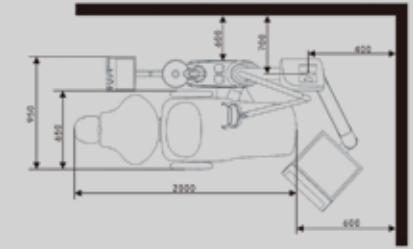

External floor-box installation sketch map

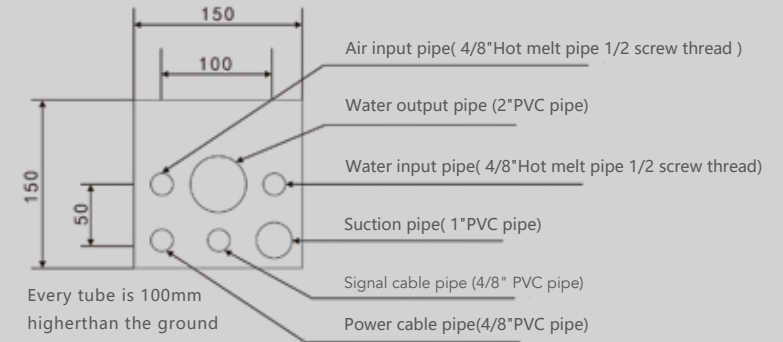

Pre-buried pipe distribution map

## Colour options

\* Please refer to actual color swatches for upholstery leather colors.



> Popular colours ★★★★★

★★★★ S06 ★★ S01 ★ S09

Fibre leather



# CINGOL

Concern about International dental Needs,  
our Goal is your Obliging Laugh.



Remark:

We reserve the right to improve all the products  
without notice in advance;  
The product's picture is only for your reference;  
the copyright is reserved and belongs to Cingol.

FOSHAN CINGOL MEDICAL INSTRUMENT CO., LTD  
Tel:0086-757-81797099 Fax:0086-757-81797055  
Email:cingol@cingol.com QQ:2851731555  
Site:www.cingol.com  
Add: Building 6-1, Nansha Industrial Village,Jiangli North Rd,  
Danzao Town,Nanhai District,Foshan,China.